

Clontech Laboratories, Inc.

Mir-X™ miRNA First-Strand Synthesis and SYBR® qRT-PCR User Manual

Cat. Nos. 638313, 638314, 638315 & 638316
(112315)

Clontech Laboratories, Inc.

A Takara Bio Company

1290 Terra Bella Avenue, Mountain View, CA 94043, USA

U.S. Technical Support: tech@clontech.com

United
States/Canada
800.662.2566

Asia Pacific
+1.650.919.7300

Europe
+33.(0)1.3904.6880

Japan
+81.(0)77.543.6116

Page 1 of 8

Table of Contents

I. Introduction.....	3
II. List of Components.....	4
III. Additional Materials Required.....	4
IV. Polyadenylation and cDNA Synthesis	5
A. Good PCR Practices.....	5
B. RNA Preparations Suitable for miRNA Analysis.....	5
C. Protocol: DNase I Treatment of RNA (Optional)	5
D. Protocol: Polyadenylation and Reverse Transcription.....	6
V. Quantification of miRNA by qPCR	6
A. Choosing Your miRNA-Specific 5'-primer for qPCR.....	6
B. Protocol: Quantifying miRNA by qPCR.....	7
VI. Data Analysis and Calculating miRNA Levels.....	7
A. Delta-Delta C _t Method	8
B. Absolute Quantification Method (Standard Curve)	8

Table of Figures

Figure 1. Mir-X miRNA qRT-PCR SYBR Kits use a single-step, single-tube reaction to produce first-strand cDNA, which is then specifically and quantitatively amplified using a miRNA-specific primer and SYBR Advantage qPCR chemistry 3

Table of Tables

Table 1. DNase I Reaction	5
Table 2. Poly(A)/cDNA Synthesis Reaction.....	6
Table 3. Sample qPCR Reaction.....	7
Table 4. U6 qPCR Reaction.....	7
Table 5. Sample Data for Delta-Delta C _t Calculation of miRNA Levels	8

Contact Us	
Customer Service/Ordering	Technical Support
tel: 800.662.2566 (toll-free)	tel: 800.662.2566 (toll-free)
fax: 800.424.1350 (toll-free)	fax: 800.424.1350 (toll-free)
web: www.clontech.com	web: www.clontech.com
e-mail: orders@clontech.com	e-mail: tech@clontech.com

I. Introduction

MicroRNAs (miRNAs) are members of a rapidly growing class of small non-coding RNAs, 19–23 nt in length, known to contribute to the regulation of gene expression in plants and animals. First discovered in *C. elegans*, miRNAs have been subsequently shown to be encoded in the genomes of higher vertebrates, including humans. miRNAs act on target mRNAs in a **sequence-specific manner** to either **promote their cleavage** and **degradation** or **reduce their translational efficiency**.

The **Mir-X miRNA qRT-PCR SYBR Kit** contains the components needed to quantify miRNAs (or other RNAs) isolated from any source. The **Mir-X miRNA First-Strand Synthesis Kit** (included in the Mir-X miRNA qRT-PCR SYBR Kit and also sold separately) is used for converting miRNAs, and other RNA molecules, into cDNA to enable specific RNAs to be quantified by real-time PCR (Figure 1). In a simple, single-tube reaction, RNA molecules are polyadenylated and reverse transcribed using poly(A) polymerase and our SMART® MMLV Reverse Transcriptase; both of which are included in the **mRQ Enzyme Mix**. The **SYBR Advantage® qPCR Premix and mRQ 3' Primer** are then used in real-time qPCR, along with your miRNA-specific 5' primer(s), to quantify specific miRNA sequences in the cDNA. Any other RNA species present in the original RNA sample, such as the mRNA targets of the miRNAs, can also be amplified and quantified using specific 5' primers.

The Mir-X miRNA kits are exceptionally easy to use and require minimal hands on time, making them perfect for high-throughput applications. Fewer manipulations means less chance of errors and highly reproducible results. The entire miRNA quantification procedure (cDNA synthesis and qPCR) can be completed within a half a day.

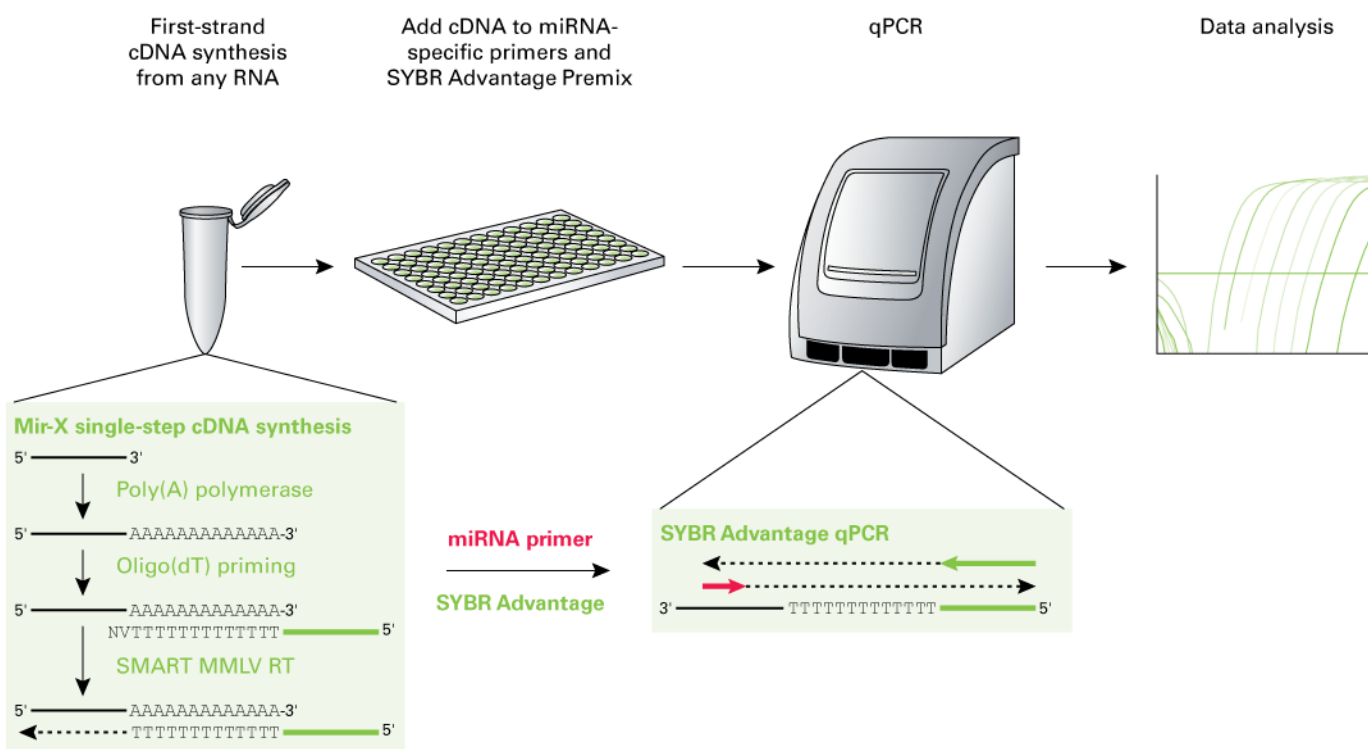


Figure 1. Mir-X miRNA qRT-PCR SYBR Kits use a single-step, single-tube reaction to produce first-strand cDNA, which is then specifically and quantitatively amplified using a miRNA-specific primer and SYBR Advantage qPCR chemistry. In the Mir-X cDNA synthesis reaction, RNAs are poly(A)-tailed using poly(A) polymerase, and then copied using a modified oligo(dT) primer and SMART MMLV Reverse Transcriptase.

II. List of Components

Store all components at -20°C .

Mir-X miRNA First-Strand Synthesis Kit (20 rxn; Cat. No. 638313)

- 25 μl mRQ Enzyme Mix
- 100 μl mRQ Buffer (2X)
- 250 μl mRQ 3' Primer (10 μM)
- 50 μl U6 Forward Primer (10 μM)
- 50 μl U6 Reverse Primer (10 μM)
- Mir-X miRNA First-Strand Synthesis and SYBR qRT-PCR User Manual (PT4445-1)

Mir-X miRNA First-Strand Synthesis Kit (60 rxn; Cat. No. 638315)

- 3 each Mir-X miRNA First-Strand Synthesis Kit (20 rxn; Cat. No. 638313)

Mir-X miRNA qRT-PCR SYBR Kit (200 rxn; Cat. No. 638314)

- 1 each Mir-X miRNA First-Strand Synthesis Kit (20 rxn; Cat. No. 638313)
 - 25 μl mRQ Enzyme Mix
 - 100 μl mRQ Buffer (2X)
 - 250 μl mRQ 3' Primer (10 μM)
 - 50 μl U6 Forward Primer (10 μM)
 - 50 μl U6 Reverse Primer (10 μM)
 - Mir-X miRNA First-Strand Synthesis and SYBR qRT-PCR User Manual (PT4445-1)
- 1 each SYBR Advantage qPCR Premix (200 rxn; Cat. No. 639676)
 - 625 μl SYBR Advantage qPCR Premix (2X; 4 each)
 - 100 μl ROX Reference Dye LSR (50X)
 - 100 μl ROX Reference Dye LMP (50X)
 - SYBR Advantage qPCR Premix User Manual (PT3883-1)

Mir-X miRNA qRT-PCR SYBR Kit (600 rxn; Cat. No. 638316)

- 3 each Mir-X miRNA First-Strand Synthesis Kit (20 rxn; Cat. No. 638313)
- 3 each SYBR Advantage qPCR Premix (200 rxn; Cat. No. 639676)

III. Additional Materials Required

- Work areas and pipettors free of contaminating DNA, RNA, and nucleases
- 0.5 M EDTA
- Cloned DNase I (RNase-free) (Takara Bio USA, Cat. No. 2270A) [Optional]
- Microfuge tubes, 0.2 ml, RNase-free
- RNase-free water
- Thermocycler, for incubating cDNA synthesis reactions
- Real-time qPCR instrument (e.g. Stratagene Mx3000P, or equivalent)

IV. Polyadenylation and cDNA Synthesis

Please read each protocol completely before starting. Successful results depend on understanding and performing all the steps correctly.

A. Good PCR Practices

Because the cDNA produced by these protocols is intended to be used in qPCR analyses, due care must be taken to ensure that the samples and reagents are not contaminated with extraneous DNA from the laboratory environment. Due to the tremendous amplification power and sensitivity of qPCR, even trace amounts of contaminating DNA and RNA will be amplified and will affect C_t and final copy number values. Before you begin, prepare work areas free of potentially contaminating DNA, RNA, and nucleases. Ideally, assemble the cDNA synthesis reactions and dilute samples and controls in one work area with a dedicated set of pipettors. Assemble the qPCR reactions in a separate area or in a noncirculating containment hood using a different set of dedicated pipettors. Wear gloves at all times and use PCR pipette tips with hydrophobic filters, and dedicated solutions. We also recommend setting up negative template control (NTC) reactions lacking any template. Finally, perform all post-PCR analyses in a separate area, preferably in a separate room, with different pipettors.

B. RNA Preparations Suitable for miRNA Analysis

Because miRNAs are very small and not polyadenylated, not all RNA preparation methods will yield miRNAs. Acceptable sources of miRNAs include: premade total RNAs available from Clontech (e.g. any of our mouse tissue total RNAs); small RNA purified using the NucleoSpin miRNA kit (Cat. Nos. 740971.10, 740971.50 & 740971.250); or total RNA purified using TRIzol or other phenol-based methods. All of these miRNA preparations can be used for miRNA cDNA synthesis and qPCR. However, spin kits or columns used for total or poly(A)+ RNA isolation are not acceptable methods of preparing miRNA, as the small size of the RNAs allow them to pass through any column.

C. Protocol: DNase I Treatment of RNA (Optional)

If your RNA sample was prepared from transfected cells, or contains trace amounts of plasmid or significant amounts of genomic DNA, the RNA should be treated with DNase I to reduce or eliminate background signal in subsequent qPCR applications. This is especially important if the miRNA of interest was expressed from a transfected plasmid or a viral vector. Contaminating DNA fragments that contain your miRNA sequence will be amplified in control reactions lacking reverse transcriptase (No RT controls), which can lead to inaccurate qPCR results. For analyses of endogenous miRNAs, DNase I treatment may not be needed.

1. Cloned DNase I (RNase-free) is available from Takara Bio USA (Cat. No. 2270A).
2. In a 0.2 ml RNase-free tube, combine the following:

Table 1. DNase I Reaction

Reagent	Volume (μl)
RNA sample (3.5–10 μg)	≤ 44
DNase I Buffer (10x)	5
DNase I (5 U/μl)	0.5
RNase-free ddH ₂ O	44.5 [RNA μl]
Total volume	50

3. In a thermocycler, incubate 30 min at 37°C, then inactivate the DNase I by adding 1 μl 0.5M EDTA and heating at 80°C for 2 min.
4. The DNase I-treated RNA can be used directly in the Polyadenylation and Reverse Transcription Protocol (Section D).

D. Protocol: Polyadenylation and Reverse Transcription

Prepare a cDNA synthesis reaction for each RNA sample you wish to analyze by qPCR. In addition, if you plan to determine the absolute level of an miRNA using a standard curve, you will need to generate cDNA from dilutions prepared from a known concentration of synthetic miRNA.

1. In an RNase-free 0.2 ml tube, combine the following reagents:

Table 2. Poly(A)/cDNA Synthesis Reaction

Reagent	Volume (μl)
mRQ Buffer (2x)	5
RNA sample (0.25–8 μg)	3.75
mRQ Enzyme	1.25
Total volume	10

2. Preferably in a thermocycler, incubate the tube for 1 hour at 37°C, then terminate at 85°C for 5 min to inactivate the enzymes.
3. Add 90 μl ddH₂O to bring the total volume to 100 μl.
4. Your cDNA is now ready for the miRNA quantification protocols in Section V.

V. Quantification of miRNA by qPCR

Before performing qPCR, determine the quantification method to be used for the samples you wish to analyze. We suggest using either a comparative C_t method (e.g. the delta-delta C_t method or ddCt), or an absolute quantification method that entails constructing a standard curve. Each of these methods requires performing additional qPCR amplifications: U6 snRNA controls for the ddC_t method, or cDNA prepared from synthetic miRNA for the standard curve method, and are described briefly in Section VI.

For additional information regarding the set-up of qPCR reactions using the SYBR Advantage qPCR Premix, including instrument-specific instructions, consult the SYBR Advantage qPCR Premix User Manual (PT3883-1).

A. Choosing Your miRNA-Specific 5'-primer for qPCR

We recommend that the entire sequence of your mature miRNA (21–23 nt) be used as your miRNA-specific, 5' primer. The 3' primer for qPCR is the mRQ 3' Primer supplied with the kit.

NOTE: The miRBase Sequence Database is a searchable database of miRNA sequences, maintained by the University of Manchester at: <http://www.mirbase.org>

B. Protocol: Quantifying miRNA by qPCR

For miRNA quantification, we recommend using the delta-delta C_t method to determine the level of each miRNA relative to the level of U6 snRNA. For this method, you will need to perform a qPCR amplification of U6 for each cDNA sample you wish to analyze. Alternatively, absolute quantification can be performed by using the C_t value(s) obtained for the miRNA in the unknown sample (or sample dilution) to obtain copy number values from a C_t vs copy number standard curve that is generated from cDNA synthesized from tenfold serial dilutions of a known concentration of synthetic miRNA.

1. Perform all sample, U6, and standard curve reactions in duplicate. Be sure to include appropriate no template controls (NTC) for each primer set. For each reaction, combine the following reagents:

Table 3. Sample qPCR Reaction

Reagent	Volume (μl)
ddH ₂ O	9
SYBR Advantage Premix (2X)	12.5
ROX Dye (50X)	0.5
miRNA-specific Primer (10 μM)	0.5
mRQ 3' Primer	0.5
cDNA	2.0
Total volume	25

Table 4. U6 qPCR Reaction

Reagent	Volume (μl)
ddH ₂ O	9
SYBR Advantage Premix (2X)	12.5
ROX Dye (50X)	0.5
U6 Forward Primer (10 μM)	0.5
U6 Reverse Primer (10 μM)	0.5
cDNA	2.0
Total volume	25

2. Cycle the reactions according to the instructions provided by the manufacturer of your real-time qPCR instrument. We routinely use the following conditions with the Stratagene Mx3000P instrument:
 - Denaturation
 - 95°C 10 sec
 - qPCR x 40 Cycles
 - 95°C 5 sec
 - 60°C 20 sec
 - Dissociation Curve
 - 95°C 60 sec
 - 55°C 30 sec
 - 95°C 30 sec

VI. Data Analysis and Calculating miRNA Levels

The two quantification approaches cited here are described only very briefly. For detailed explanations of these and other qPCR strategies, consult a comprehensive qPCR reference. For additional information regarding the set-up of qPCR reactions using the SYBR Advantage qPCR Premix, including instrument-specific instructions, consult the SYBR Advantage qPCR Premix User Manual (PT3883-1).

A. Delta-Delta C_t Method

The delta-delta C_t method (ddCt) is an approximation method that measures the relative levels of a miRNA between two samples by comparing them to a second RNA which serves as a normalization standard (e.g. U6). Briefly, the unknown miRNA and the U6 RNA are amplified in each sample to determine the C_t for each. This allows relative levels to be determined using the ddCt calculation. An example is shown below:

Table 5. Sample Data for Delta-Delta C_t Calculation of miRNA Levels

Sample RNA	Ct(miR-122a)	Ct(U6)	dCt
Brain	30	25	5
Liver	20	24	-4

Calculating relative copy numbers using the delta-delta C_t method:

$$[\text{miR-122a}]_{\text{Liver}}/[\text{miR-122a}]_{\text{Brain}}: 2^{-\text{dCt}(\text{Liver})} / 2^{-\text{dCt}(\text{Brain})} = 2^{-(-4)} / 2^{-(5)} = 29 = 512^*$$

*This result indicates that the level of miR-122a is 512-fold higher in the liver than in the brain.

Theoretically, between successive PCR cycles (i.e., each delta C_t) there is a twofold difference in PCR product. In this case, a normalized C_t difference of 9 between the samples gives 2⁹, which is 512-fold. **Important:** a number of assumptions need to be considered before this method can be used correctly. The relative levels of the U6 should not differ between samples. More importantly, the amplification efficiency of the primers should be near 100% (i.e., the amount of product should actually increase twofold with each PCR cycle). If primer efficiency is significantly less than 100%, then this method will not be accurate and either the primers should be redesigned or the absolute quantification method should be used (Section B).

B. Absolute Quantification Method (Standard Curve)

In this method, a calibrated preparation of synthetic miRNA is used to make serial dilutions and generate the standard curve. The plot then used to determine miRNA copy number from the C_t values obtained from the experiment samples.

1. Determine the average C_t values for the duplicate qPCR reactions of the cDNA samples generated from the diluted synthetic miRNA samples, and plot these C_t values vs. input miRNA copy number on a log scale.
2. Determine average C_t values for each duplicate experiment sample, or sample dilution, and extrapolate the corresponding RNA copy number from these C_t values using the standard curve generated in Step 1. Use all C_t values that are below that of the NTC.

Notice to Purchaser

Our products are to be used for research purposes only. They may not be used for any other purpose, including, but not limited to, use in drugs, *in vitro* diagnostic purposes, therapeutics, or in humans. Our products may not be transferred to third parties, resold, modified for resale, or used to manufacture commercial products or to provide a service to third parties without prior written approval of Clontech Laboratories, Inc.

Your use of this product is also subject to compliance with any applicable licensing requirements described on the product's web page at <http://www.clontech.com>. It is your responsibility to review, understand and adhere to any restrictions imposed by such statements.

SYBR is a registered trademark of Molecular Probes, Inc.

Clontech®, the Clontech logo, Advantage, Mir-X, and SMART are trademarks of Clontech Laboratories, Inc. All other marks are the property of their respective owners. Certain trademarks may not be registered in all jurisdictions. Clontech is a Takara Bio Company. ©2015 Clontech Laboratories, Inc.

This document has been reviewed and approved by the Clontech Quality Assurance Department.